


OVERVIEW: ICX connection screws




ICX hexagonal instrument
SW 1.2mm

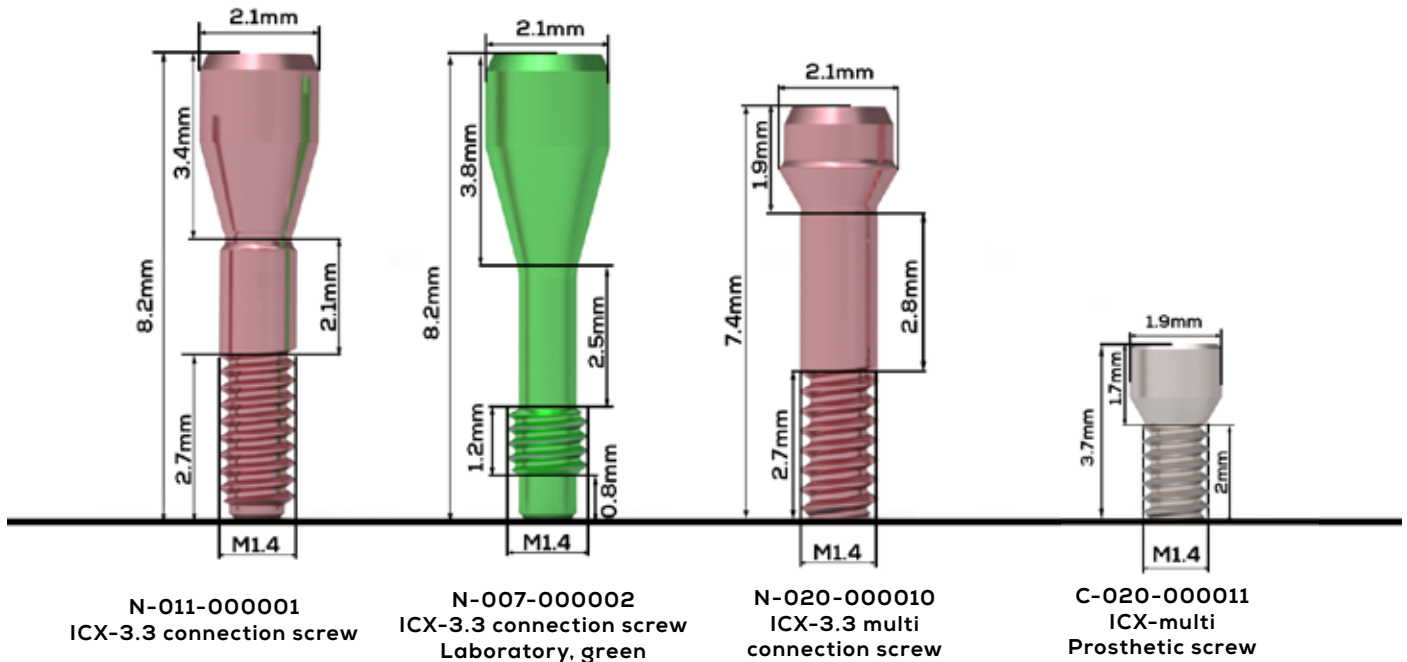


OVERVIEW: ICX connection screws



ICX hexagonal instrument
SW 1.4mm

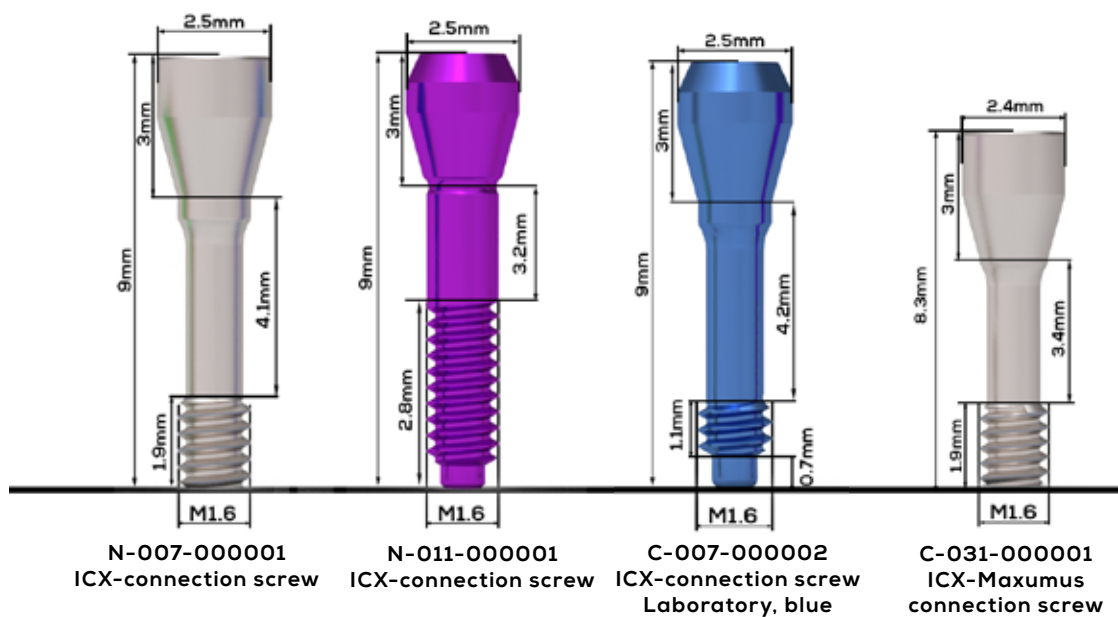

ICX hexagonal instrument
SW 1.2mm

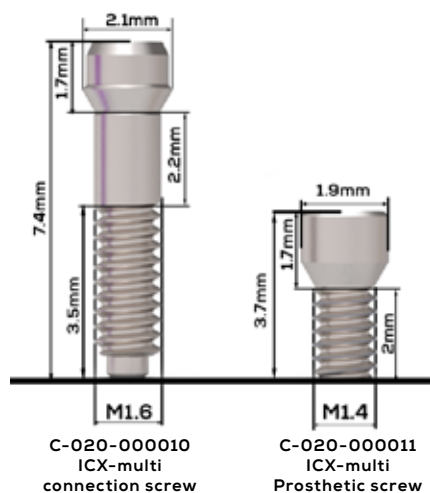

RECOMMENDED VALUES FOR TORQUES



5 Ncm to 10 Ncm (hand-tight)	
Screw plugs	
Gingiva former	
Impression posts	
Plastic abutments	

RECOMMENDED VALUES FOR TORQUES



max. 15 Ncm	
Scan body and all other abutments made of PEEK	A collection of various PEEK abutments and scan bodies. The components are shown in different shapes and sizes, including single abutments, stacked abutments, and scan bodies with different base shapes.

RECOMMENDED VALUES FOR TORQUES



max. 30 Ncm	
all titanium abutments	
Adhesive bases	
CERICX	
Ball head	
Maximus	
ICX-multi abutments	