

ICX

THE FAIR
PREMIUM
IMPLANT-
SYSTEM.

We're not number one – but we're on our way there!

Combined experience, taken a step further – with a clear commitment to quality and precision.


ICX-FLEX ADHESIVE BASES:

Maximum flexibility for the angled screw channel.

The ICX-Flex adhesive base with an angled screw channel offers an innovative solution for challenging cases involving an unfavourable implant position or aesthetically demanding conditions.




0°-30° angulation
Can be shortened diagonally & vertically!

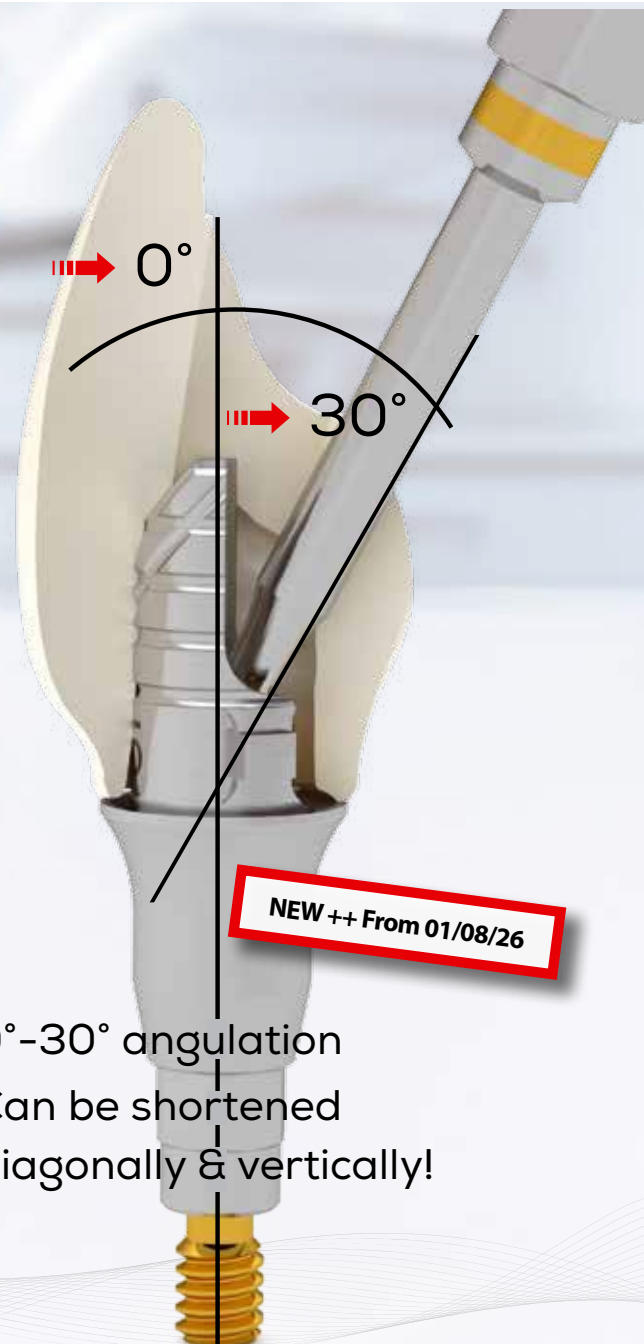

ONE VISION – BE SMART. BE ICX.

medentis[®]
medical

THE ICX-FLEX TITANIUM ADHESIVE BASES.

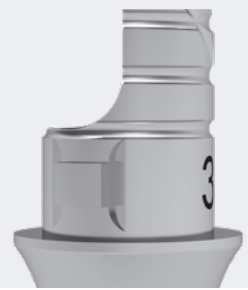


THE ICX-FLEX ADHESIVE BASE WITH AN ANGLED SCREW CHANNEL – PRECISION AND FLEXIBILITY FOR THE HIGHEST STANDARDS

The ICX-Flex adhesive base combines innovative technology and the highest quality to offer dentists and dental technicians a reliable and adaptable solution.

Thanks to its sophisticated design, the ICX-FLEX adhesive base allows for precise adjustment and alignment of the screw channel. This can be flexibly positioned within the mouth and individually shortened to two different lengths. This ensures optimal functionality and aesthetics, even for particularly demanding requirements, especially in the anterior region.


ICX-FLEX GINGIVAL HEIGHTS

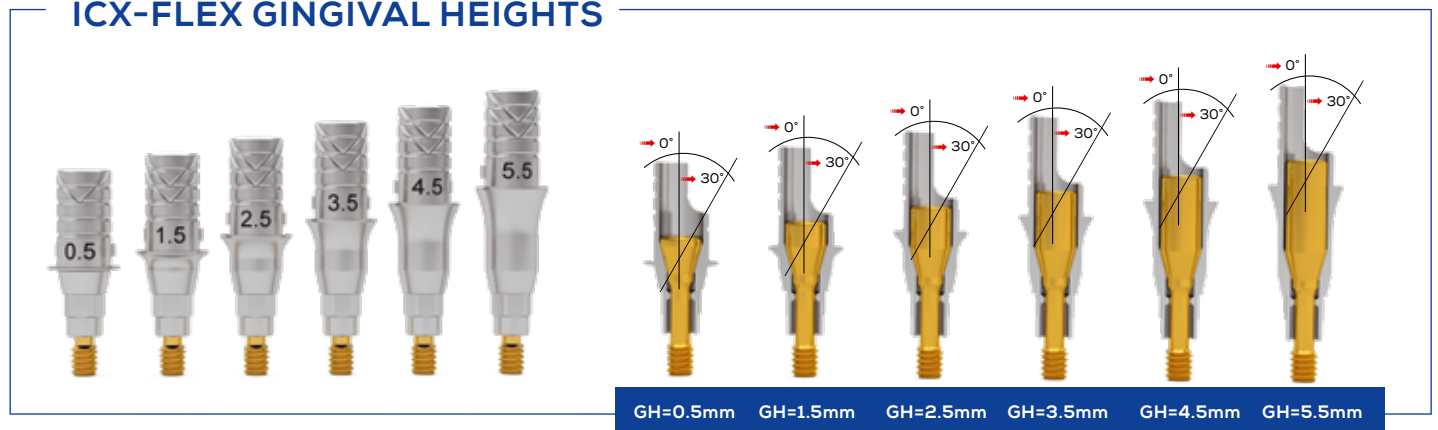

IMPORTANT NOTE: All screws for the Flex adhesive bases have a Torx head and are **ALWAYS** gold anodized. Furthermore, each screw corresponds to a specific gingival height and is not interchangeable with others.

ADVANTAGES OF THE ICX-FLEX ADHESIVE BASES

- can be shortened diagonally & vertically
- angled screw channel
- double anti-rotation protection
- optimised adhesive interface
- optimised and slimmer emergence profile

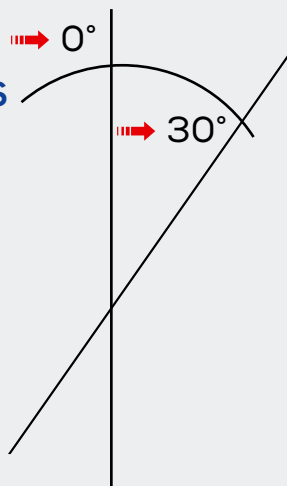

→ 30° angle



medentis
medical

→ Can be shortened diagonally & vertically!

THE ICX-FLEX TITANIUM ADHESIVE BASES.

VARIABLE CHIMNEY HEIGHT (CH)

The abutment height (AH) of 6 mm supports even tall restorations; however, it can be individually shortened diagonally to 5 and 4 mm, thereby adapting perfectly to the clinical situation where vertical clearance is limited.

OPTIMISED EMERGENCE PROFILE

The optimised emergence profile supports and protects the soft tissue.

ANGLED SCREW CHANNEL

The ball-Torx screwdriver socket allows for angles of up to 30°.


Horizontal
CH=6

Can be shortened diagonally CH=5
Can be shortened diagonally CH=4

0°

30°

NEW ++ From 01/08/26

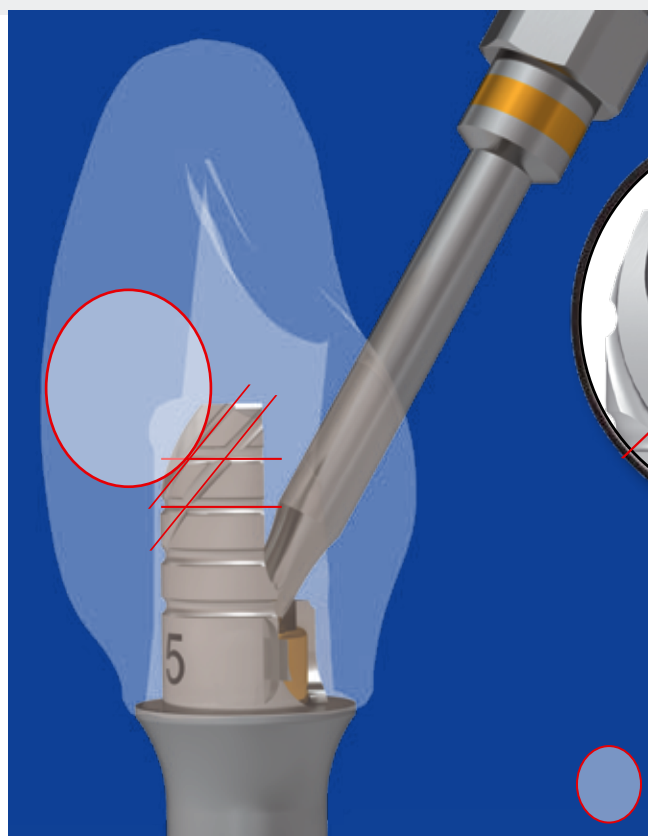

CH=5

CH=4

0°-30° angle

Can be shortened diagonally!
For optimal wall thickness
of the zirconia crowns.




ICX-Flex in our shop!

FOR IMPROVED PROSTHETIC RESTORATIONS


The ICX-Flex adhesive base for ICX-Standard and ICX-3.3mm expands the existing adhesive base portfolio with additional functional features. It features a straight screw channel and can be used either as a rotating or non-rotating component, thereby supporting various restoration options for single crowns and bridges.

The design of the ICX-Flex adhesive base incorporates various gingival heights. The longer chimney is designed so that it can be individually adjusted to the desired height if required.

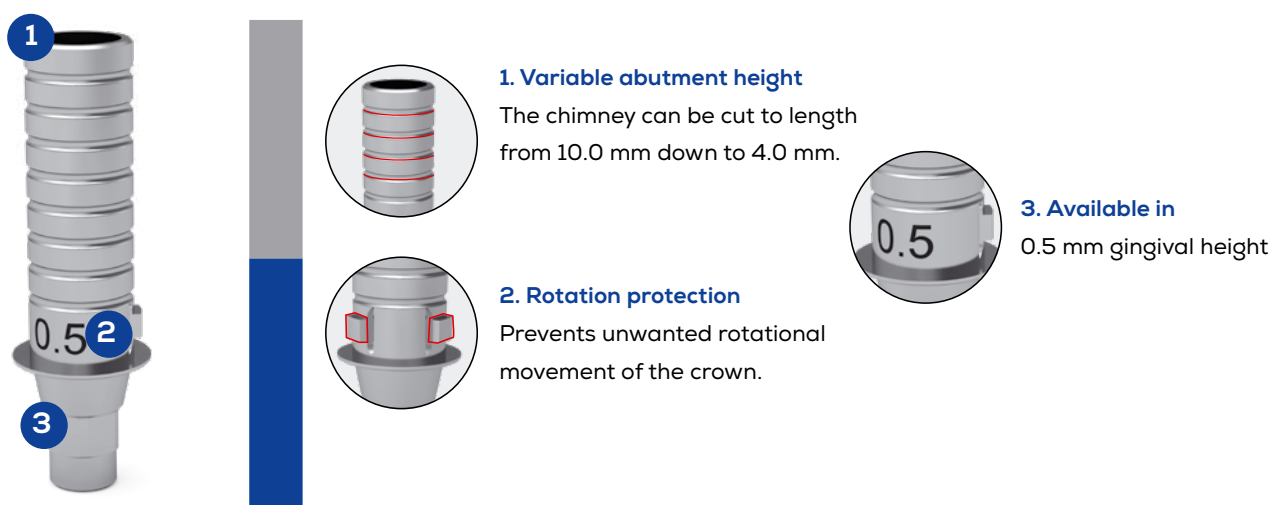

NON-HEX ADHESIVE BASE FOR BLOCKED SUPRA-STRUCTURES

Our non-hex adhesive abutments with a high collar have been specially developed for fixed implant-supported superstructures. By eliminating the anti-rotational hex connection, they allow for a less stressful and more passive fit – ideal for multi-unit bridges and complex implant axes.

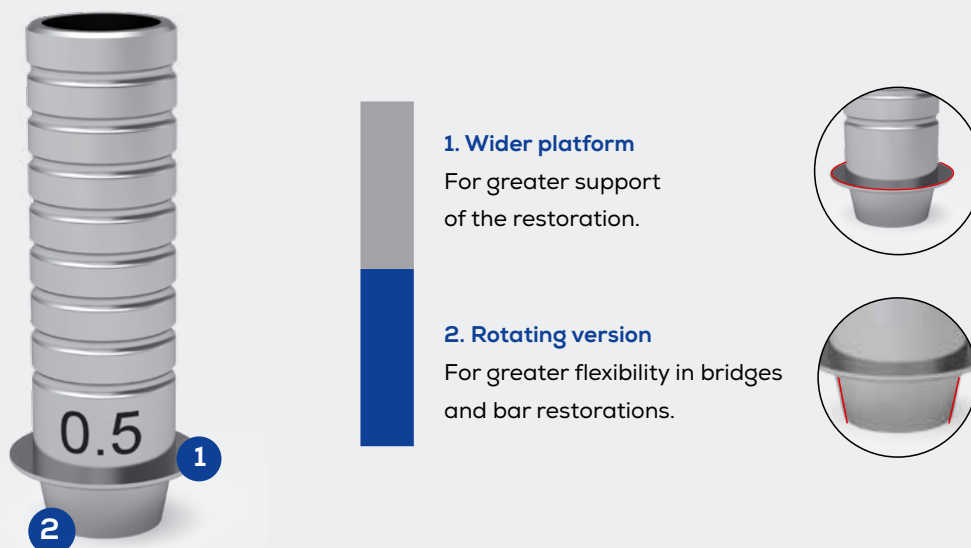


THE ICX-FLEX TITANIUM ADHESIVE BASES.

FOR CASES WITH HIGH AESTHETIC DEMANDS

The innovative design of the ICX-Flex adhesive base has been specially developed for angled screw channels and is ideally suited to the requirements of prosthetic restorations in the aesthetic zone.


ICX-FLEX ADHESIVE BASE STRAIGHT AND ANGLED

In the event of an unfavourable implant position or an implant position in the aesthetic zone, it is possible to relocate the screw channel to a non-visible area.

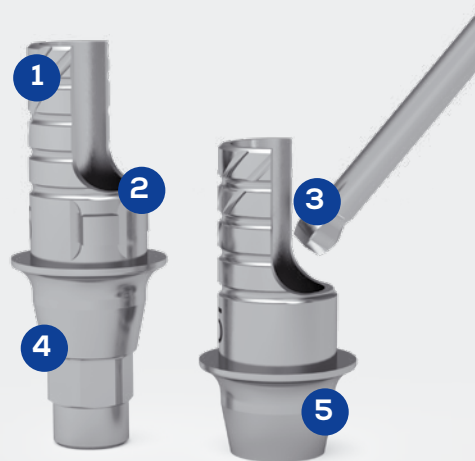

1. Adjustable chimney height

The chimney can be individually shortened from 6.0 mm to 5.0 mm and 4.0 mm.

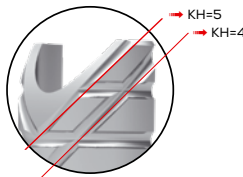

2. Screw channel

Compatible with angled and straight screw channels.



3. Ball-Torx insertion instrument

For angulations up to 30°.
For effortless and reliable screw tightening.



4. Available in several gingival heights

Gingival heights ranging from 0.5 mm to 5.5 mm for greater flexibility.



5. Rotating version

For greater flexibility with bridges and bar restorations.



NEW ++ From 01/08/26



ICX-Flex in our shop!

IMPORTANT NOTE::

**THE APPLICATIONS AND PRODUCT FEATURES ARE IDENTICAL
FOR THE ICX-FLEX STANDARD AND 3.3MM ADHESIVE BASES!**

ICX-FLEX TITANIUM ADHESIVE BASES

OVERVIEW: ICX-FLEX System Ø3.3mm



ICX-Flex in our shop!

NEW ++ From 01/08/26

ICX-Flex adhesive bases 3.3 mm

Article-no.	Article 3.3mm	Dimensions/material		Price in €
N-029-011101	ICX-Flex 3.3 adhesive base	GH=0.5mm	hex	45,-
N-029-011107	ICX-Flex 3.3 adhesive base	GH=1.5mm		45,-
N-029-011108	ICX-Flex 3.3 adhesive base	GH=2.5mm		45,-
N-029-011109	ICX-Flex 3.3 adhesive base	GH=3.5mm		45,-
N-029-011104	ICX-Flex 3.3 adhesive base	GH=4.5mm		45,-
N-029-011105	ICX-Flex 3.3 adhesive base	GH=5.5mm		45,-
N-029-011001	ICX-Flex 3.3 adhesive base, non-hex	GH=0.5mm	non hex	45,-
N-029-011007	ICX-Flex 3.3 adhesive base, non-hex	GH=1.5mm		45,-
N-029-011008	ICX-Flex 3.3 adhesive base, non-hex	GH=2.5mm		45,-
N-029-011009	ICX-Flex 3.3 adhesive base, non-hex	GH=3.5mm		45,-
N-029-011004	ICX-Flex 3.3 adhesive base, non-hex	GH=4.5mm		45,-
N-029-011005	ICX-Flex 3.3 adhesive base, non-hex	GH=5.5mm		45,-
N-029-201101	ICX-Flex, 0°, 3.3 adhesive base, hex	K = 10mm, GH 0,5mm	hex non hex	45,-
N-029-201001	ICX-Flex, 0°, 3.3 adhesive base, non-hex	K = 10mm, GH 0,5mm		45,-
C-015-100050	Screwdriver, Torx	short		49,-
C-015-100051	Screwdriver, Torx	medium		49,-
C-015-100052	Screwdriver, Torx	long		49,-
N-030-000006	ICX Scan Base 2.0	New!		39,-

ICX-Flex 3.3mm Adhesive Bases / Healing Caps

Article-no.	Article 3.3mm	Dimensions/material		Price in €
N-004-362001	ICX-Flex 3.3 healing cap, straight	Ø3.6 K=2 for GH=0.5mm		29,-
N-004-362007	ICX-Flex 3.3 healing cap, straight	Ø3.6 K=2 for GH=1.5mm		29,-
N-004-362008	ICX-Flex 3.3 healing cap, straight	Ø3.6 K=2 for GH=2.5mm		29,-
N-004-362009	ICX-Flex 3.3 healing cap, straight	Ø3.6 K=2 for GH=3.5mm		29,-
N-004-362004	ICX-Flex 3.3 healing cap, straight	Ø3.6 K=2 for GH=4.5mm		29,-
N-004-362005	ICX-Flex 3.3 healing cap, straight	Ø3.6 K=2 for GH=5.5mm		29,-




ICX-Flex in our shop!

OVERVIEW:

ICX-FLEX System - Standard Ø3.45 · Ø3.75 · Ø4.1 · Ø4.8mm

NEW ++ From 01/08/26

ICX-Flex adhesive bases standard Ø3.45 · 3.75 · 4.1 · 4.8mm

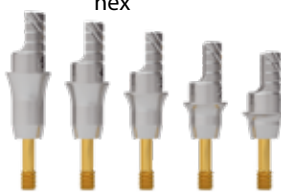Article-no.	Article Standard Ø3.45 · 3.75 · 4.1 · 4.8mm	Dimensions/material		Price in €
C-029-011102	ICX-Flex adhesive base, standard	GH=0.5mm		45,-
C-029-011107	ICX-Flex adhesive base, standard	GH=1.5mm		45,-
C-029-011108	ICX-Flex adhesive base, standard	GH=2.5mm		45,-
C-029-011109	ICX-Flex adhesive base, standard	GH=3.5mm		45,-
C-029-011104	ICX-Flex adhesive base, standard	GH=4.5mm		45,-
C-029-011105	ICX-Flex adhesive base, standard	GH=5.5mm		45,-
C-029-011002	ICX-Flex adhesive base, standard non-hex	GH=0.5mm		45,-
C-029-011007	ICX-Flex adhesive base, standard non-hex	GH=1.5mm		45,-
C-029-011008	ICX-Flex adhesive base, standard non-hex	GH=2.5mm		45,-
C-029-011009	ICX-Flex adhesive base, standard non-hex	GH=3.5mm		45,-
C-029-011004	ICX-Flex adhesive base, standard non-hex	GH=4.5mm		45,-
C-029-011005	ICX-Flex adhesive base, standard non-hex	GH=5.5mm		45,-
C-029-201102	ICX-Flex adhesive base, standard	GH=0.5mm; K=10mm		45,-
C-029-201002	ICX-Flex adhesive base, standard non-hex	GH=0.5mm; K=10mm		45,-
C-015-100050	Screwdriver, Torx	short		49,-
C-015-100051	Screwdriver, Torx	medium		49,-
C-015-100052	Screwdriver, Torx	long		49,-
C-030-000006	ICX Scan Base 2.0	New!		39,-

ICX-Flex Standard Adhesive Bases / Healing Caps

Article-no.	Article Standard Ø3.45 · 3.75 · 4.1 · 4.8mm	Dimensions/material		Price in €
C-004-442002	ICX-Flex healing cap, straight, standard	Ø4.4 K=2 for GH=0.5mm		29,-
C-004-442007	ICX-Flex healing cap, straight, standard	Ø4.4 K=2 for GH=1.5mm		29,-
C-004-442008	ICX-Flex healing cap, straight, standard	Ø4.4 K=2 for GH=2.5mm		29,-
C-004-442009	ICX-Flex healing cap, straight, standard	Ø4.4 K=2 for GH=3.5mm		29,-
C-004-442004	ICX-Flex healing cap, straight, standard	Ø4.4 K=2 for GH=4.5mm		29,-
C-004-442005	ICX-Flex healing cap, straight, standard	Ø4.4 K=2 for GH=5.5mm		29,-
C-004-652002	ICX-Flex healing cap, conical, standard	Ø6.5 K=2 for GH=0.5mm		29,-
C-004-652007	ICX-Flex healing cap, conical, standard	Ø6.5 K=2 for GH=1.5mm		29,-
C-004-652008	ICX-Flex healing cap, conical, standard	Ø6.5 K=2 for GH=2.5mm		29,-
C-004-652009	ICX-Flex healing cap, conical, standard	Ø6.5 K=2 for GH=3.5mm		29,-
C-004-652004	ICX-Flex healing cap, conical, standard	Ø6.5 K=2 for GH=4.5mm		29,-
C-004-652005	ICX-Flex healing cap, conical, standard	Ø6.5 K=2 for GH=5.5mm		29,-

ICX

THE FAIR
PREMIUM
IMPLANT-
SYSTEM.

Proven expertise, applied efficiently.

ICX-FLEX ADHESIVE BASES:


Maximum flexibility for the angled screw channel.

Proven expertise,
consistently applied –
for predictable and
reliable restorations.

NEW ++ From 01/08/26



ICX-Flex in our Shop!

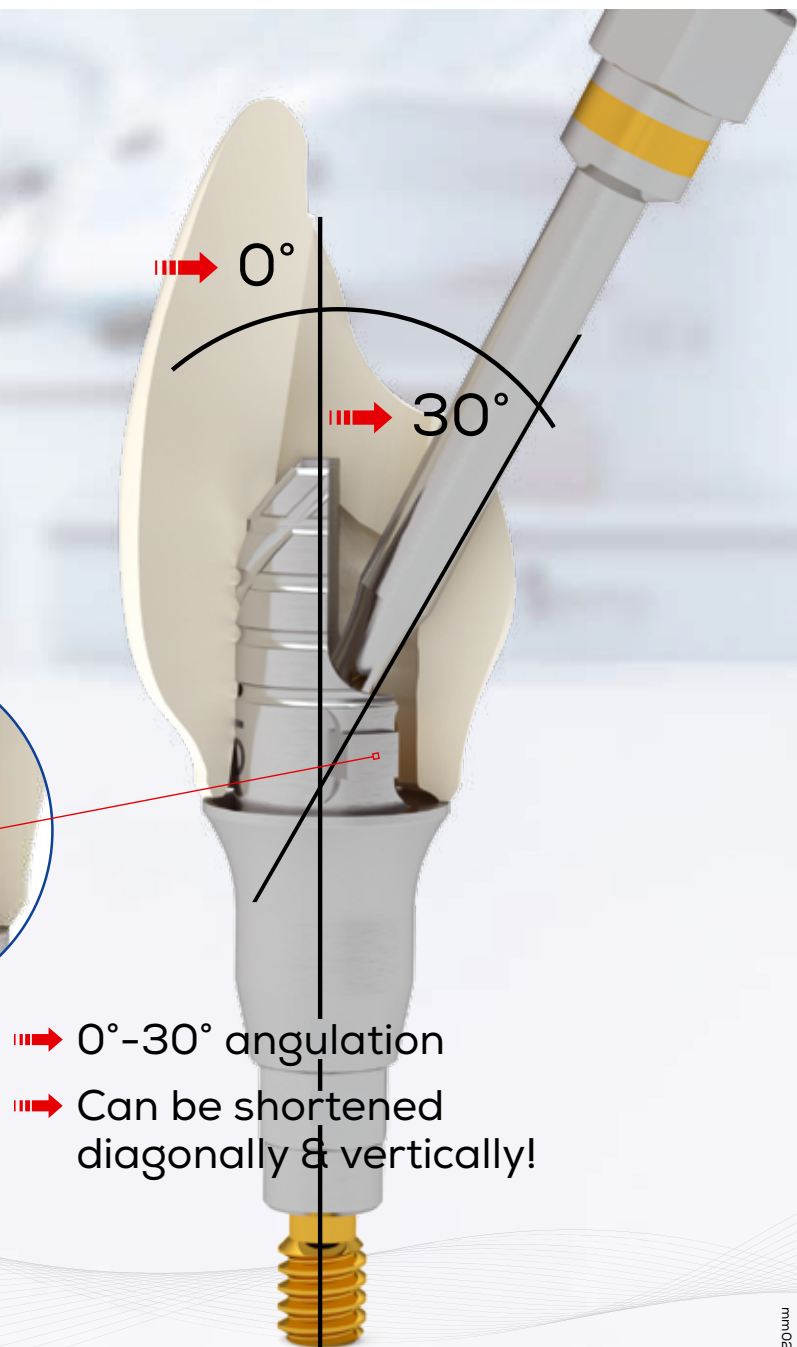

- 0°-30° angulation
- Can be shortened diagonally & vertically!

mm02846

ONE VISION – BE SMART. BE ICX.

medentis[®]
medical